Datasheet

HCA-S-200M-IN

200 MHz Photoreceiver with InGaAs PIN Photodiode



The picture shows the HCA-S-200M-IN-FS with free space input. The photoreceiver will be delivered without post holder and post.

Features	 InGaAs PIN detector Spectral range 900 1700 nm Bandwidth DC 200 MHz Amplifier transimpedance (gain) 2.0 x 10⁴ V/A Max. conversion gain 1.9 x 10⁴ V/W @ 1550 nm 		
Applications	 Spectroscopy Fast pulse and transient measurements Optical triggering Optical front-end for oscilloscopes, A/D converters and HF lock-in amplifiers 		
Specifications	Test conditions	$V_{_S} = \pm 15$ V, $T_{_A} = 25$ °C, system impedance = 50 Ω	
Gain	Amplifier transimpedance Max. conversion gain	2.0 x 10 ⁴ V/A (@ 50 Ω load) 1.9 x 10 ⁴ V/W (@ 1550 nm)	
Frequency Response	Lower cut-off frequency Upper cut-off frequency (–3 dB) Rise/fall time (10 % - 90 %)	DC 200 MHz (±15 %) 1.8 ns	
Detector	Detector material Active area Spectral response	InGaAs PIN photodiode \varnothing 300 μ m (free space "-FS" version only) 900 1700 nm	
Input	Input offset compensation range Optical saturation power NEP	$\pm 100~\mu A$ adjustable by offset potentiometer 60 μW (for linear amplification, @ 1550 nm) 5.2 pW/ \sqrt{Hz} (@ 1550 nm, 10 MHz)	

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E M T O

Datasheet

Specifications (continued)

HCA-S-200M-IN

200 MHz Photoreceiver with InGaAs PIN Photodiode

Output Output voltage range $\pm 1.2 \, \text{V}$ (@ 50 Ω load)

for linear amplification and low harmonic distortion

Max. output voltage range $\pm 1.7 \text{ V}$ (@ 50 Ω load)

Output impedance 50Ω (designed for 50Ω load)

Output noise typ. 30 mV_{PP} or 4.5 mV_{RMS}

(@ 50 Ω load, no signal on detector)

Power Supply Supply voltage $\pm 15 \text{ V}$

Supply current ± 60 mA typ.

(depends on operating conditions, recommended

power supply capability min. ±150 mA)

Case Weight 210 g (0.5 lbs)


Material AlMg4.5Mn, nickel-plated

Temperature Range Storage temperature -40 ... +100 °C

Operating temperature 0 ... +60 °C

Absolute Maximum Ratings Optical input power 10 mW Power supply voltage ±22 V

Spectral Response

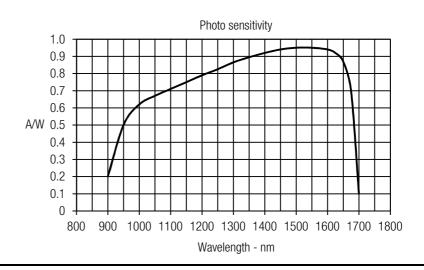

Connectors

Input HCA-S-200M-IN-FS 25 mm round flange for

free space applications


HCA-S-200M-IN-FC FC fiber optic receptacle

Output BNC jack (female)

Power supply Lemo® series 1S, 3-pin fixed socket

(Mating plug type: FFA.1S.303.CLAC52)

Pin 1: +15 V Pin 2: -15 V Pin 3: GND

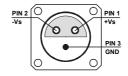


SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E MT O

Typical Performance Characteristics Frequency response



Noise spectrum

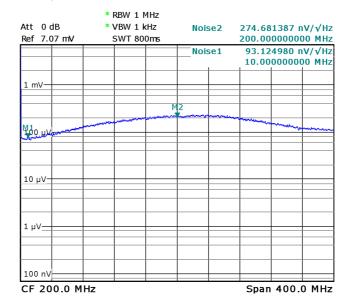


Note: Spectral noise data is measured at the amplifier output with no signal on the photodiode. To determine the spectral input noise divide the measured output noise by the amplifier conversion gain.


Conversion gain (V/W) = amplifier gain (20,000 V/A) x photo sensitivity (A/W).

Marker	Frequency	Output noise	Resulting input noise (NEP)
1	10 MHz	93 nV/√Hz	4.9 pW/√Hz (@ 1550 nm)
2	200 MHz	275 nV/√Hz	15 pW/√Hz (@ 1550 nm)

Typical Performance Characteristics (continued) Pulse response to square wave input signal (with 16 times averaging)

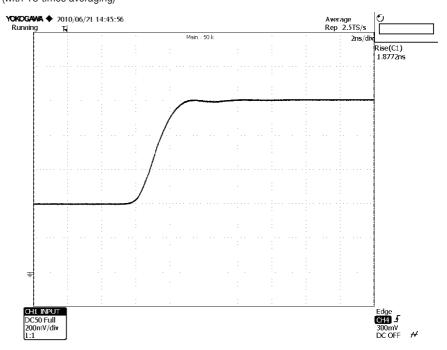


Large signal response output signal for 100 MHz, 55 μW modulated optical input signal (with 4 times averaging)

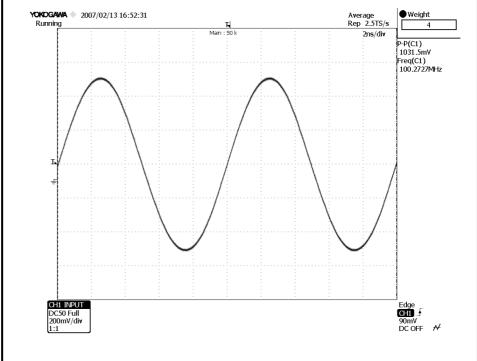


SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E M T O

Typical Performance Characteristics (continued) Small signal response output signal for 1.2 μ W modulated optical input signal, 1 MHz square wave (without (top) and with 64 times averaging (bottom))

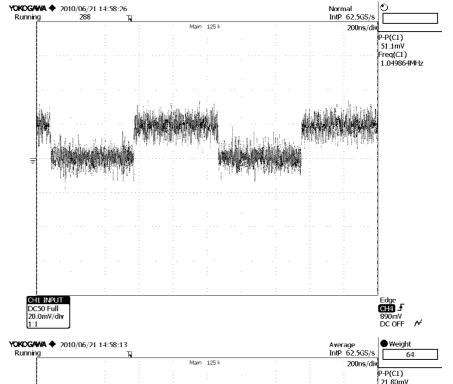

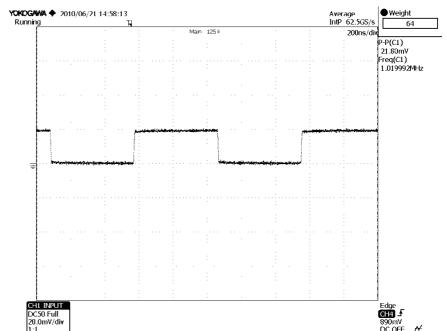

Available Models


HCA-S-200M-IN-FS HCA-S-200M-IN-FC HCA-S

free space input FC fiber optic receptacle customized versions available on request


SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

Е М ТО

Dimensions

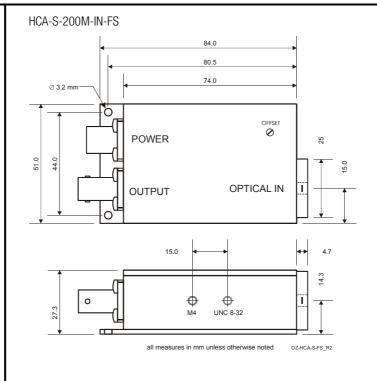

HCA-S-200M-IN-FC

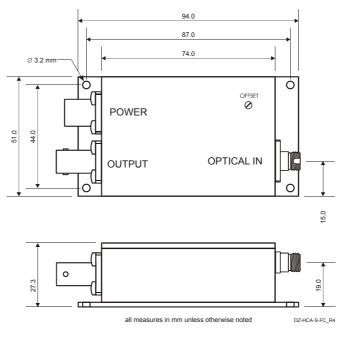

FEMTO Messtechnik GmbH Klosterstr. 64 10179 Berlin · Germany Phone: +49 30 280 4711-0 Fax: +49 30 280 4711-11 Email: info@femto.de www.femto.de Specifications are subject to change without notice. Information provided herein is believed to be accurate and reliable. However, no responsibility is assumed by FEMTO Messtechnik GmbH for its use, nor for any infringement of patents or other rights of third parties which may result from its use. No license is granted by implication or otherwise under any patent or patent rights of FEMTO Messtechnik GmbH. Product names mentioned may also be trademarks used here for identification purposes only.

© by FEMTO Messtechnik GmbH · Printed in Germany

F E T O