

TUNABLE LASER SOURCE

UP TO 20W INPUT POWER



The tunable laser source (TLS) provides continuous output from 400 nm to 1000 nm (VIS) and 1000 nm to 2500 nm (SWIR) with bandwidths (FWHM) of 2.5 nm and 5 nm respectively. Custom and extended spectral ranges (up to 2500 nm) and bandwidths (sub-nm) are also available. Photon etc. tunable Laser Source (TLS) is compatible with any VIS-NIR broadband source. This high-end product provides the highest out-of-band reflection (< -60 db) available on the market. It is an ideal tool for instruments calibration, spectroscopy and hyperspectral imaging.

TECHNICAL SPECIFICATIONS

	CONTRAST VIS	CONTRAST SWIR	CONTRAST X
Spectral Range (extended and reduced spectral ranges also available**)	400-1000 nm	1000-2300 nm (2500 nm optional*)	where X represents a custom spectral range
Bandwidth (FWHM)***	1.0 - 2.5 nm	2.0 - 5.0 nm	High resolution 0.15 nm - 0.9 nm
Out of Band Rejection (fibered output or with the background suppressor accessory)	< -60 dB @ ± 40 nm	< -60 dB @ ± 80 nm (measured up to 1.7 μm)	Depends on the bandwidth
Maximum input average power	HP8 (up to 8W), HP20 (up to 20W)	HP8 (up to 8W), HP20 (up to 20W)	HP4 (up to 4W)
Peak Efficiency	typically around 65%		
Optical Density (OD)	> OD6 (measured at 1064 nm)		TBD
Damage Threshold	< 5 GW/cm2 peak power @ 1064 nm, 8 ns		
Beam Diameter	2 to 5 mm (depending on λ)		
Input Beam Divergence Requirement	< 1 mrad		
Wavelength Resolution (Relative)	FWHM / 8		
Pointing Stability	< 1 mm lateral displacement @ 1 m from filter		
Scanning speed (multiple step)	35 ms stabilization time for 0.1 nm step, 45 ms stabilization time for 0.2 nm step, 55 ms stabilization time for 1 nm step, 60 ms stabilization time for 2 nm step, 65 ms stabilization time for 5 nm step, 70 ms stabilization time for 10 nm step		
Operating System (OS)	Windows Vista (32 & 64 bits), Windows 7 (32 & 64 bits), Windows 8 (32 & 64 bits)		
Software	PhySpec™ included (SDK available)		
Computer Connection	USB 2.0 (compatible 1.1)		
Dimensions (L x W x H)	9" x 6.3" x 6.7" / 23 cm x 16 cm x 17 cm (filter)		
Operating Temperature	5 to 40°C		
Storage Temperature	0 to 50°C		
Power Supply	100 - 240 V , 50 - 60 Hz		


OPTIONS & ACCESSORIES

Enhance SWIR	N/A	* up to 2500 nm	
Fibered Output	An X-Y-Z translation adjustment allows coupling optimization		
Harmonic Filter	Blocks the harmonics coming from the region 400-500 nm	Blocks the harmonics coming from the region 500-1000 nm and/or 1000-1250 nm	Filter chosen according to spectral range
Alignment Kit (for free space)	In free-space configuration, the alignment kit allows the user to rapidly find the correct alignment		
** Extended and reduced spectral ranges also available	e.g.: 500-2000 nm, 400-1700nm, 500-900 nm, 400-650 nm, 650-1000 nm, 1000-1700 nm, 1700-2300 nm, etc.		

*** Valid if the divergence of the input beam does not exceed 1 mrad

NOTE: Photon etc reserves the right to change the design and specification of the product at any time, without notice.

LLTF VIS AND SWIR OUTPUT POWER

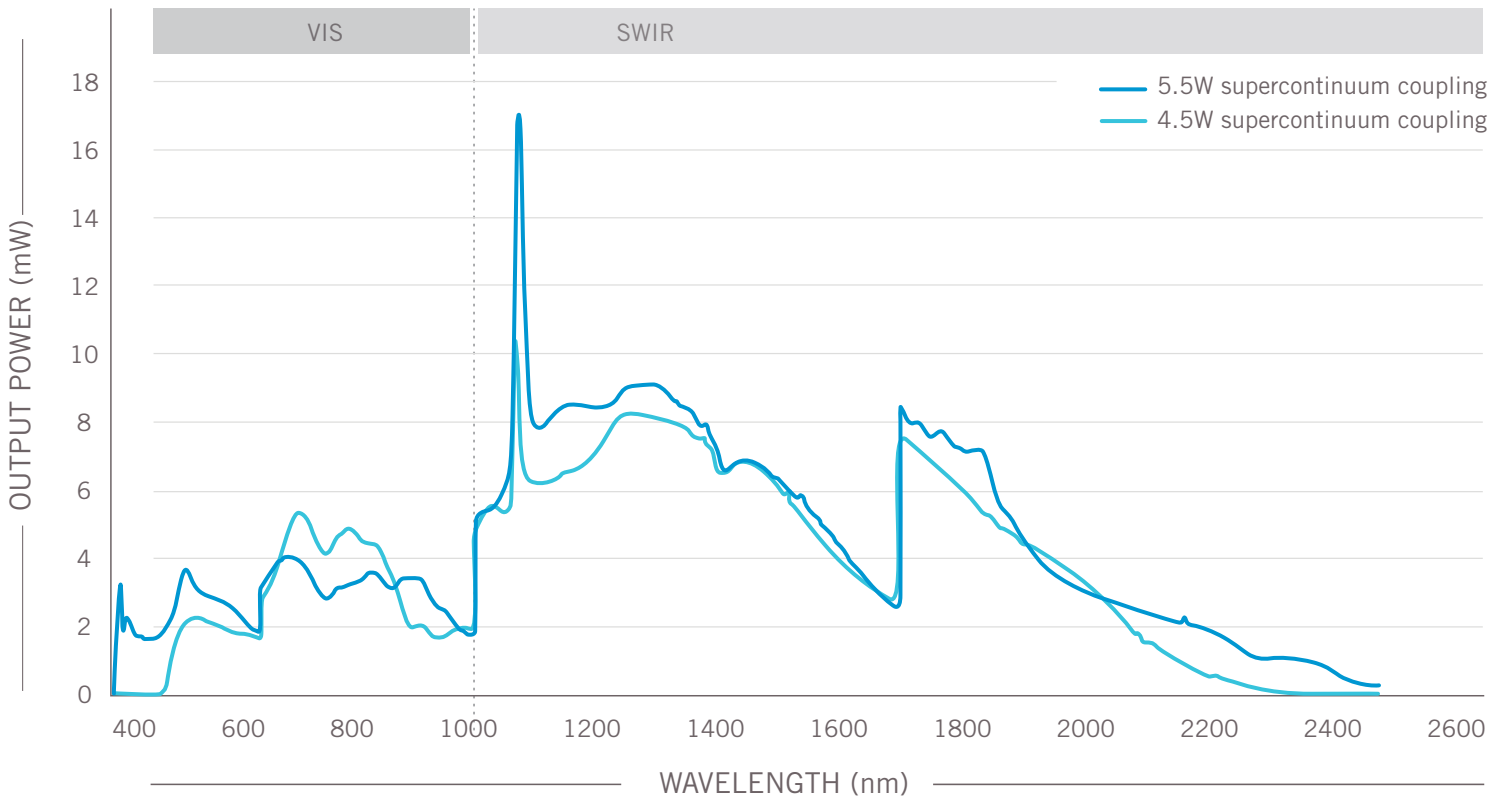

LLTF EXTENDED (500-2000 nm) OUTPUT POWER


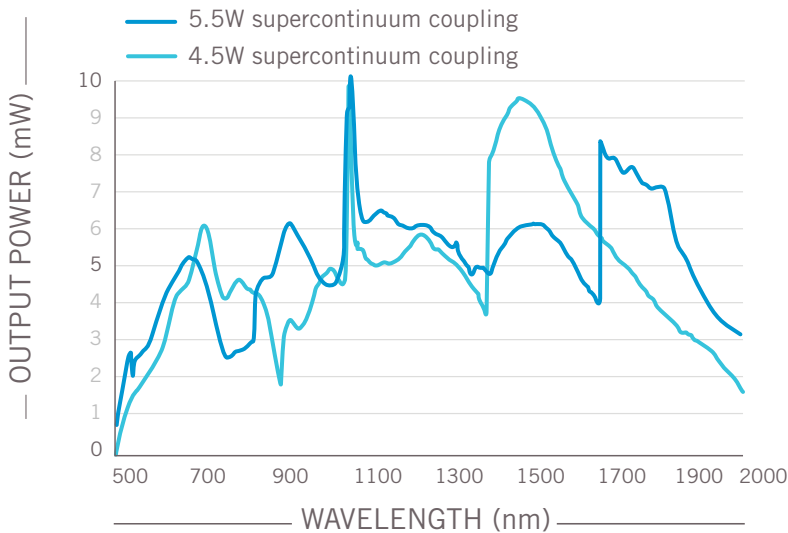


Illustration of the out-of-band rejection of a volume holographic grating at $\lambda_c = 632$ nm. Bands of ± 45 nm are presented and an out-of-band rejection of -60 dB is obtained.

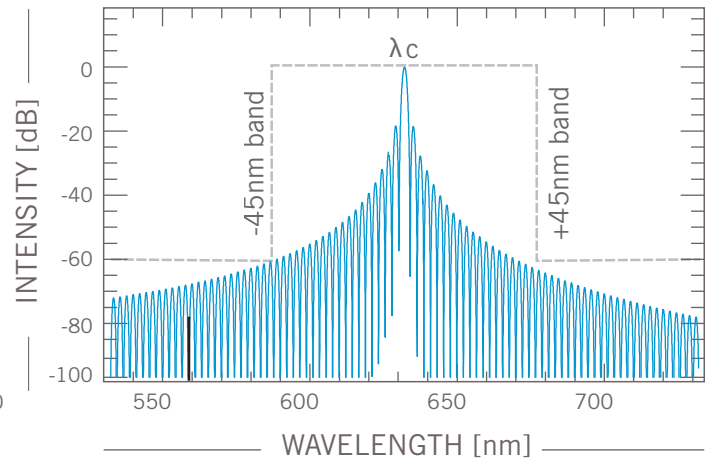

OUR DISTRIBUTORS

EUROPE

ITALY
CRISEL INSTRUMENTS
+39 06 35402933
vitucci@crisel-instruments.it

PORTUGAL
IZASA SCIENTIFIC LDA
+351 21 424 73 00
info-pt@izasascientific.com

BELGIUM, NETHERLANDS & LUXEMBURG
TE LINTELO SYSTEMS B.V.
+31 316-340804
roland@tlsbv.nl

UNITED KINGDOM & IRELAND
PRO-LITE TECHNOLOGY LTD
+44 (0) 12344 36110
nick.barnett@pro-lite.co.uk

GERMANY, SWITZERLAND, AUSTRIA
MD INNOVATION TECH GMBH
+49 9233 7157745
info@md-innovationtech.com

SPAIN
IZASA SCIENTIFIC S.L.U.
+902 20 30 80
jrecasens@izasascientific.com

FRANCE
OPTON LASER INTERNATIONAL
+33 (0)1 69 41 04 05
ventes@optonlaser.com

POLAND
LABNATEK
+22 119 47 17
info@labnatek.pl

GERMANY, SWITZERLAND, AUSTRIA & EASTERN EUROPE
SPHEREOPTICS
+49 (0)8152 983 78-90
info@sphereoptics.de

SWEDEN, NORWAY, DENMARK & FINLAND
KMG FOTONIK AB
+46 (0) 70 64 64 275
info@kmgfotonik.se

ASIA & OCEANIA

JAPAN
TOKYO INSTRUMENT
+81 3 3686 4711
sales@tokyoinst.co.jp

SINGAPOUR
CADENCE TECHNOLOGIES
+65 6779 0900
cwchong@cadence-tech.com.sg

AUSTRALIA
SCITECH PTY LTD
+61 3 9480 4999
con@scitech.com.au

CHINA
AUNION TECH
+86 21 51083793
info@auniontech.com

SOUTHEAST ASIA
ALV TECHNOLOGIES
+65 67468068

INDIA
INKARP INSTRUMENTS PVT.LTD
+914027172293
harishk@inkarp.co.in

SHANGAI CHEER INSTRUMENT LTD
18201727656

GENERAL TECH SCIENCES
+86 10 68337300
info@gtechsci.com