Datasheet FWPR-20-SI

Femtowatt Photoreceiver with Si Photodiode



The photoreceiver will be delivered without post holder and post.

Features	 Si photodiode, 1.1 x 1.1 mm² active area Ultra low noise, NEP 0.7 fW/√Hz Amplifier transimpedance gain 1 x 10¹² V/A Max. conversion gain 0.6 x 10¹² V/W @ 960 nm Wavelength range 320 1100 nm Fluorescence measurements Spectroscopy Electrophoresis Replacement for photomultiplier tubes (PMTs) and avalanche photodiodes (APDs) 		
Applications			
Specifications	Test conditions	$V_{\rm S}=\pm 15$ V, $T_{\rm A}=25^{\circ}{\rm C}$ Warm-up 20 minutes (min. 10 minutes recommended)	
Gain	Amplifier transimpedance Max. conversion gain	1.0 x 10 ¹² V/A 0.6 x 10 ¹² V/W	(@ ≥ 100 k Ω load) (@ 960 nm)
Frequency Response	Lower cut-off frequency Upper cut-off frequency (–3 dB) Rise/fall time (10 % - 90 %)	DC 20 Hz 18 ms	(±20 %) (±20 %)
Detector	Detector material Active area Spectral response	Si photodiode 1.1 x 1.1 mm ² 320 1100 nm	1
Input	Optical saturation power NEP	18 pW (for linear amplification, @ 960 nm) 0.7 fW/√Hz (@ 960 nm, 1 Hz)	

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E T O

Datasheet FWPR-20-SI

Femtowatt Photoreceiver with Si Photodiode

Specifications (continued)		
Output	Output voltage range Output impedance Offset voltage Max. output current Output noise	$\begin{array}{lll} \pm 10 \ V & (@ \geq 100 \ k\Omega \ load) \\ 50 \ \Omega & (designed \ for \geq 100 \ k\Omega \ load) \\ 0 \ V, \ adjustable \ by \ offset \ potentiometer \ within \ \pm 1.6 \ V \\ \pm 25 \ mA \\ ca. \ 40 \ mV_{PP} \ or \ 6 \ mV_{RMS} \\ (@ \geq 100 \ k\Omega \ load, \ no \ signal \ on \ detector) \end{array}$
Power Supply	Supply voltage Supply current	$\pm 15~\text{V}$ $\pm 15~\text{mA}$ typ. (depends on operating conditions, recommended power supply capability min. $\pm 50~\text{mA})$
Case	Weight Material	190 g (0.42 lbs) AIMg3/4.5Mn, nickel-plated
Temperature Range	Storage temperature $-40 \dots +100 ^{\circ}\text{C}$ Operating temperature $0 \dots +60 ^{\circ}\text{C}$	
Absolute Maximum Ratings	Optical input power Power supply voltage	10 mW ±22 V
Spectral Response	0.7 0.6 0.5 0.4 A/W 0.3 0.2 0.1 0 200 300 400	Photo sensitivity

Datasheet FWPR-20-SI

Femtowatt Photoreceiver with Si Photodiode

Connectors Input 25 mm round flange for free space applications


(fiber optic input available as customized unit)

BNC jack (female) Output

Lemo® series 1S, 3-pin fixed socket **Power Supply**

(Mating plug type: FFA.1S.303.CLAC52)


Pin 1: +15V Pin 2: -15VPin 3: GND

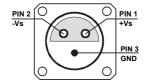

Available Models FWPR-20-SI-FS Free space input FWPR-S

Customized version available on request

Dimensions

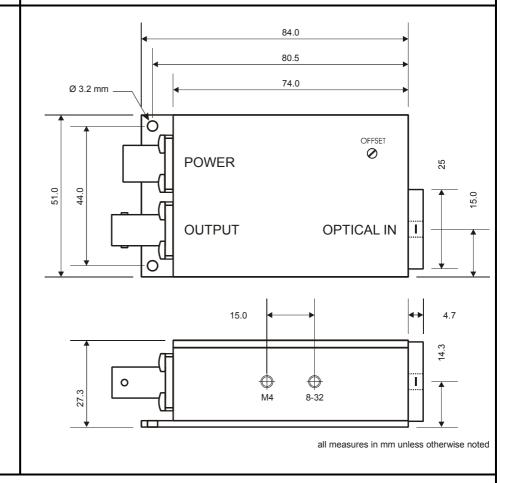

FEMTO Messtechnik GmbH Klosterstr. 64 10179 Berlin · Germany Phone: +49 30 280 4711-0 Fax: +49 30 280 4711-11 Email: info@femto.de www.femto.de

Specifications are subject to change without notice. Information provided herein is believed to be accurate and reliable. However, no responsibility is assumed by FEMTO Messtechnik GmbH for its use, nor for any infringement of patents or other rights of third parties which may result from its use. No license is granted by implication or otherwise under any patent or patent rights of FEMTO Messtechnik GmbH. Product names mentioned may also be trademarks used here for identification purposes only.

© by FEMTO Messtechnik GmbH · Printed in Germany

0